




COVIDIEN

positive results for life™

Directional Stapling Technology

Improved reliability in staple formation

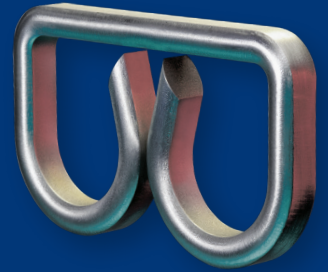


DST Series™ EEA™ Staplers and Orvil™ Devices



DST Series™
Technology

In Challenging Applications the Difference is DIRECTIONAL STAPLING TECHNOLOGY

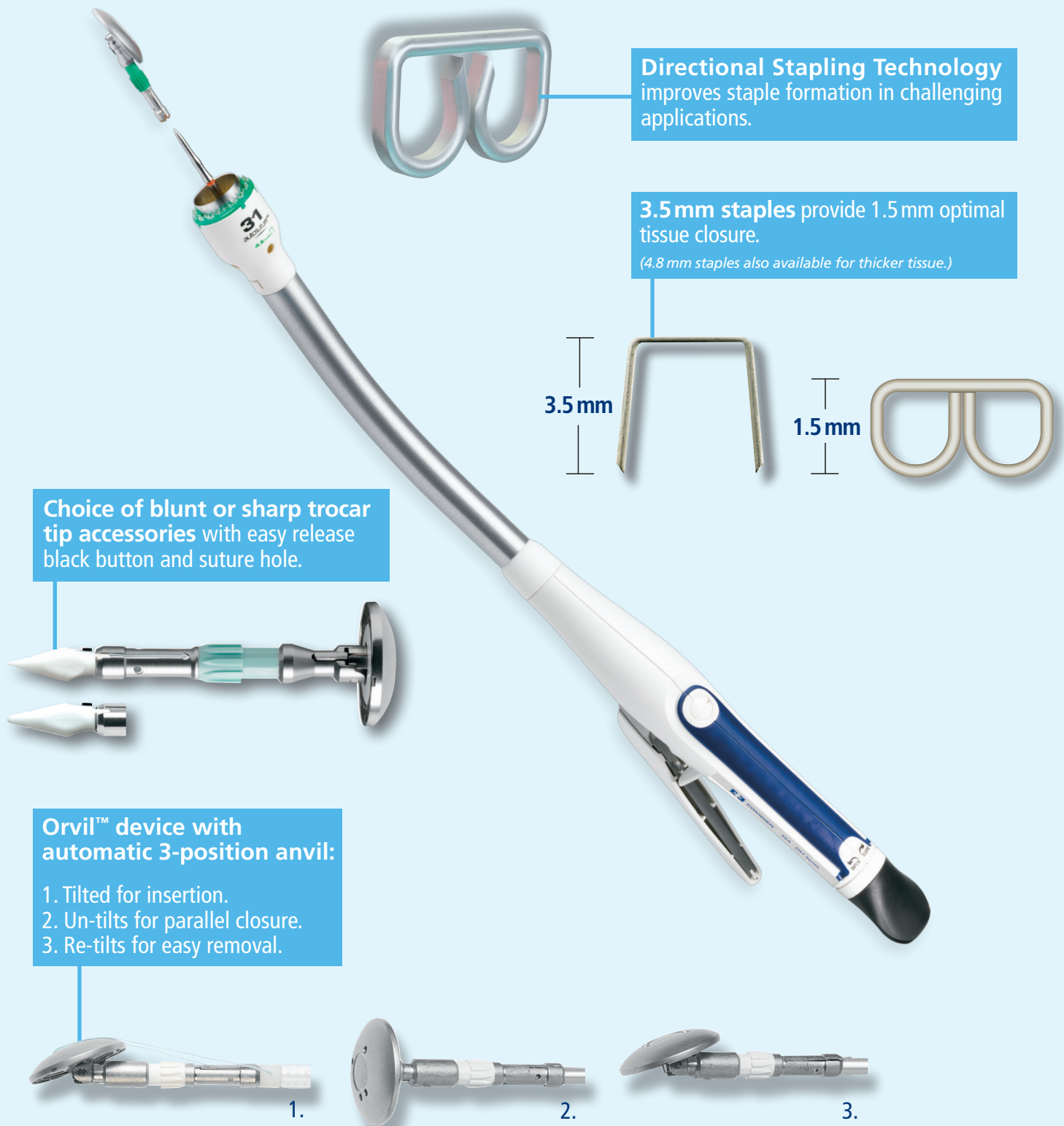

Directional Stapling Technology improves staple formation in challenging applications.

3.5 mm staples provide 1.5 mm optimal tissue closure.
(4.8 mm staples also available for thicker tissue.)

Choice of blunt or sharp trocar tip accessories with easy release black button and suture hole.

Orvil™ device with automatic 3-position anvil:

1. Tilted for insertion.
2. Un-tilts for parallel closure.
3. Re-tilts for easy removal.

Complete instructions, warnings, precautions and contraindications for this product are included in its package insert, and should be read prior to use.

In Challenging Applications the Difference is DIRECTIONAL STAPLING TECHNOLOGY



The only 25 mm stapler with 3.5 mm staples
Optimal lumen. Optimal closure.

35 cm working length
17% longer than competitive models for improved access.


Audible and tactile feedback
provides confidence at critical steps in the procedure.

First ready-to-use trans-oral delivery system delivers the anvil — fast — smooth — easy.

100 cm 18 FR delivery tube
precise delivery through a gastrotomy.





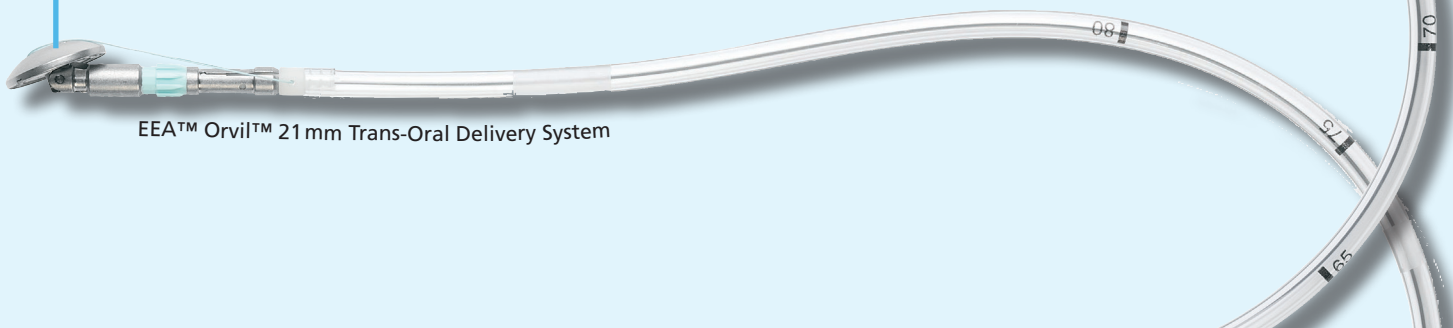





EEA™ Orvil™ 25 mm Trans-Oral Delivery System

Unique pre-tilted anvils
provide smallest cross-section for transesophageal passage.



EEA™ Orvil™ 21 mm Trans-Oral Delivery System

Stapling Made Smarter

PRODUCT OFFERING	REORDER CODE	PRODUCT DESCRIPTION	COLOR	STAPLE SIZE
	EEA21 EEAXL21	DST Series™ EEA™ 21 mm single use stapler DST Series™ EEA™ XL 21 mm single use stapler	● Aqua ● Aqua	4.8 4.8
	EEA2135 EEAXL2135	DST Series™ EEA™ 21 mm single use stapler w/3.5mm staples DST Series™ EEA™ XL 21 mm single use stapler w/3.5mm staples	● Aqua ● Aqua	3.5 3.5
	EEA25 EEAXL25	DST Series™ EEA™ 25 mm single use stapler DST Series™ EEA™ XL 25 mm single use stapler	● White ● White	4.8 4.8
	EEA2535 EEAXL2535	DST Series™ EEA™ 25 mm single use stapler w/3.5mm staples DST Series™ EEA™ XL 25 mm single use stapler w/3.5mm staples	● White ● White	3.5 3.5
	EEA28 EEAXL28	DST Series™ EEA™ 28 mm single use stapler DST Series™ EEA™ XL 28 mm single use stapler	● Blue ● Blue	4.8 4.8
	EEA2835 EEAXL2835	DST Series™ EEA™ 28 mm single use stapler w/3.5mm staples DST Series™ EEA™ XL 28 mm single use stapler w/3.5mm staples	● Blue ● Blue	3.5 3.5
	EEA31 EEAXL31	DST Series™ EEA™ 31 mm single use stapler DST Series™ EEA™ XL 31 mm single use stapler	● Green ● Green	4.8 4.8
	EEA33 EEAXL33	DST Series™ EEA™ 33 mm single use stapler DST Series™ EEA™ XL 33 mm single use stapler	● Yellow ● Yellow	4.8 4.8
	EEAORVIL21*	DST Series™ EEA™ Orvil™ 21 mm	● Aqua	N/A
	EEAORVIL25*	DST Series™ EEA™ Orvil™ 25 mm	● White	N/A

*EEAORVIL21 & EEAORVIL25 compatible with XL Staplers only.

For more information contact your local sales representative or visit www.covidien.com/surgical



COVIDIEN, COVIDIEN with logo, Covidien logo and *positive results for life* are U.S. and/or internationally registered trademarks of Covidien AG. All other brands are trademarks of a Covidien company.
© 2013 Covidien. 4.13 M120089-1

For more information or to set up a product demonstration contact your local Covidien surgical devices specialist or call customer service

555 LONG WHARF DRIVE
NEW HAVEN, CT
06511

800-722-8772

WWW.COVIDIEN.COM/SURGICAL